

 **Zymfentra**TM
infliximab-dyyb | subcutaneous
120 mg/mL

PREFILLED PEN

HOW TO USE
ZYMFENTRATM




Please read the [Medication Guide](#) for key information on serious side effects.
For more information please see complete [Instructions for Use](#).

HOW TO USE ZYMFENTRA

BEFORE YOU BEGIN

Use the prefilled pen only if you have received appropriate training from your healthcare provider and make sure you have a detailed plan for how often you will need to give an injection. Do not reuse the pen.

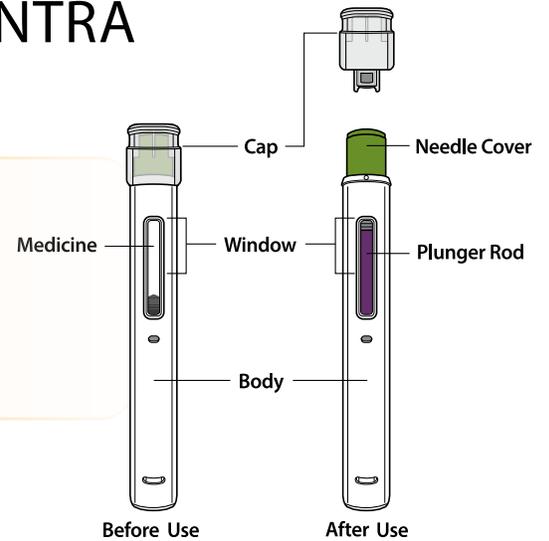


PREPARE FOR THE INJECTION

1. GATHER YOUR SUPPLIES


- Find a clean, flat, well-lit surface such as a table or countertop
- Remove the carton containing the prefilled pen from the refrigerator. If needed, the pen may also be stored at room temperature at 68°F to 77°F (20°C to 25°C) for up to 14 days
- Remove the pen from the carton
- Items needed:
 - ZYMFENTRA[™] Prefilled Pen
 - Alcohol swab
 - Cotton ball or gauze*
 - Adhesive bandage*
 - Sharps disposal container/household plastic container*

*Not included in carton.

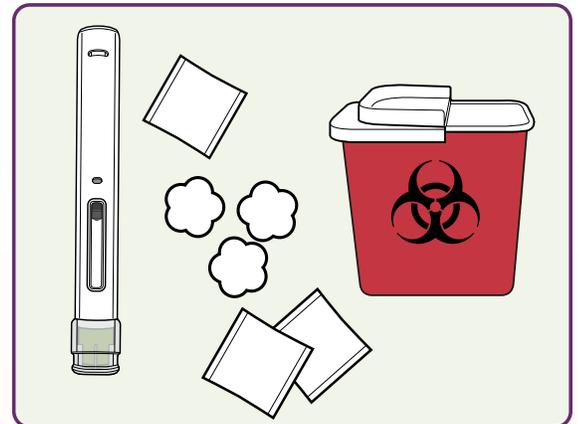

2. INSPECT THE PEN AND MEDICINE

- DO NOT use the pen if:
 - It is cracked or visibly damaged
 - The expiration date has passed
 - The liquid is discolored or contains particles (the liquid should be clear and colorless to a pale brown)
 - It has been stored at room temperature for longer than 14 days
- Do not shake the pen at any time
- It is normal to see air bubbles inside the liquid

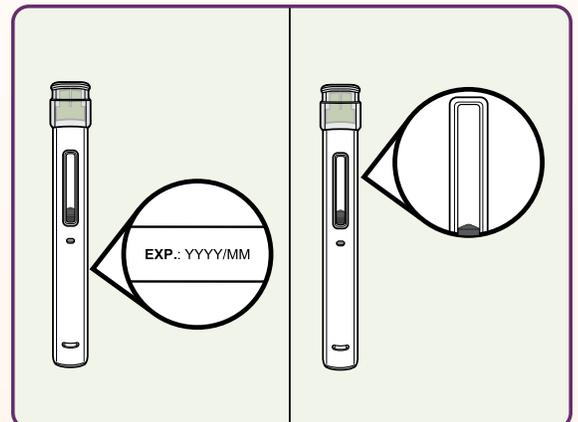


Please read the [Medication Guide](#) for key information on serious side effects. For more information please see complete [Instructions for Use](#).

PREPARE FOR THE INJECTION continued


3. WAIT 30 MINUTES

- If refrigerated, leave the pen at room temperature for 30 minutes to allow it to warm up naturally
- Do not use a microwave, stovetop, or other heat sources to warm the pen

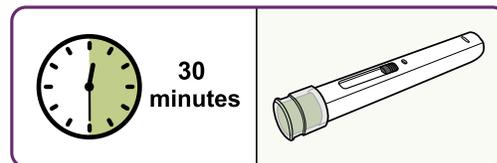


4. CHOOSE AN INJECTION SITE

- Recommended injection sites:
 - The front of the thighs
 - The abdomen
 - The outer area of the upper arm, ONLY if given by a caregiver
- Rotate the injection site each time you give an injection
- Each new injection site should be an inch or more away from the previous site
- Do not inject into the 2-inch area surrounding your belly button (navel)
- Do not inject into skin that is tender, damaged, bruised, or scarred

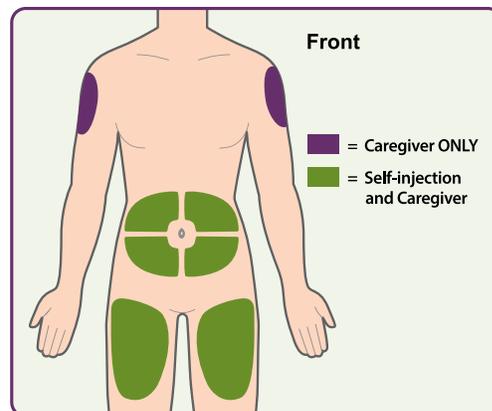

5. WASH YOUR HANDS AND CLEAN THE INJECTION SITE


- Wash your hands with soap and water and dry them fully
- Clean the injection site using the alcohol swab included in your carton
- Let your skin dry before injecting
- Do not blow on or touch the injection site again before giving the injection

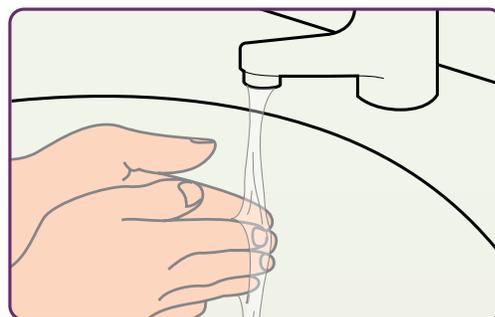

GIVE THE INJECTION


6. REMOVE THE CAP

- Remove the cap from the pen by pulling the cap straight off
- Be careful not to touch the needle cover
- Throw away the cap in your plastic container

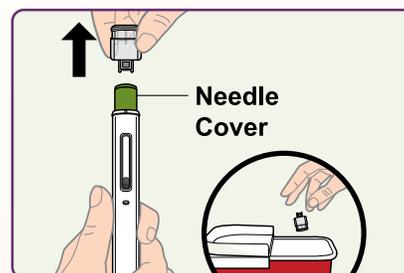

7. PLACE THE PEN

- Hold the pen so that you can see the window
- Without pinching or stretching the skin, place the pen over the injection site at a 90-degree angle

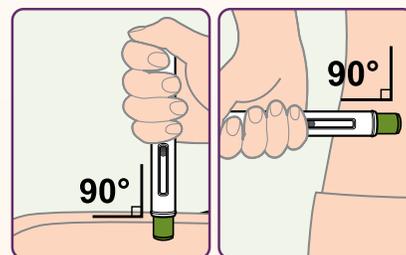

Please read the [Medication Guide](#) for key information on serious side effects. For more information please see complete [Instructions for Use](#).

 **Zymfentra**[™]
infliximab-dyyb | subcutaneous
120 mg/mL


8. START THE INJECTION

- Press the pen firmly against the skin
- When the injection starts, you will hear the first of two loud “click” sounds. The plunger rod will begin to fill the window
- Keep holding the pen and listen for the second “click” sound
- After you hear the second click, continue to hold the pen and count slowly for at least 5 seconds



9. REMOVE THE PEN


- Check that the window is filled with the plunger completely
- Lift the pen from the injection site
- If the plunger rod is not filling the window fully, you did not receive your full dose. Do not reuse the pen in this case. Call your healthcare provider immediately
- Gently press a cotton ball, gauze, or adhesive bandage over the injection site. Be sure not to rub the injection site

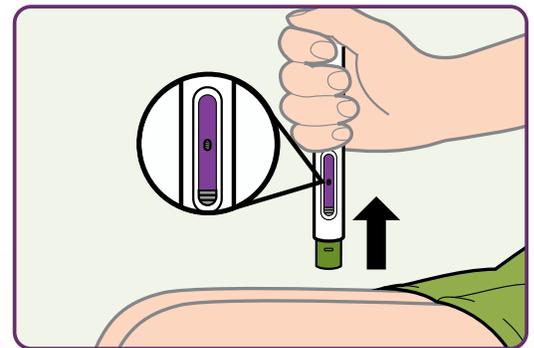

AFTER THE INJECTION

10. DISPOSE OF THE PEN

- Put the used pen in an approved sharps disposal container immediately after use. Do not recap the pen
- If you do not have an approved sharps disposal container, you may use a household container that is:
 - Made of heavy-duty plastic
 - Able to close with a tight-fitting, puncture resistant-lid
 - Properly labelled to warn of hazardous waste inside the container
- When your container is almost full, dispose of it in accordance with local requirements. For specific information about sharps disposal in your area, visit: [fda.gov/safesharpsdisposal](https://www.fda.gov/safesharpsdisposal)
- Keep the pen and sharps disposal container away from children

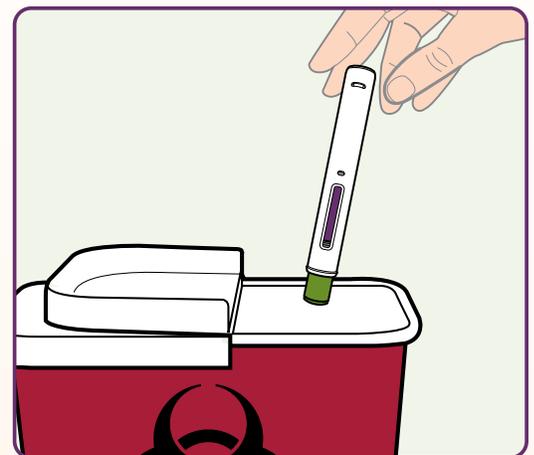

Please read the [Medication Guide](#) for key information on serious side effects. For more information please see complete [Instructions for Use](#).